COLPHARMA®

MediKilo

Bilancia pesapersona
analogica

10371

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EN

Operating Instruction Manual


 **Read the instructions carefully before using this device.**

WARNING

- This device must be used exclusively as described in this manual. This scale is designed for personal weighing and complies with international standards set for such products. The manufacturer/distributor shall not be liable for any damages caused by improper use.

HOW TO USE THE DEVICE

- Place the scale on a flat, horizontal surface. A hard surface is most suitable to obtain maximum weighing accuracy. Placing the scale on a soft surface may reduce weighing accuracy.
- Adjust the knob under the scale plate until the hand points to 0 kg.
- Step on and off the scale to make sure that the calibration is correct. If the indicator hand does not

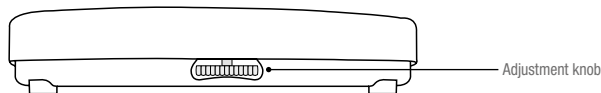return to 0 kg, repeat step 2.

- Step on the scale and stand still. Your feet should be positioned so that the weight is evenly distributed. Read the weighing value.

 **It is usual for the weight to change during the day or from one day to the other. Small differences in readings are normal.**

Cleaning/Maintenance

Use a cloth lightly moistened with a mild detergent. Please do not use thinners or benzene. No lubrication of the internal mechanisms is required.

The device must be disposed of according to local regulations and not as ordinary household waste.

TECHNICAL SPECIFICATIONS

Maximum load capacity: 160 Kg
Dimensions: 42,8 x 29 x 8 cm
Weight: 3 kg
Accuracy: 0-50 Kg \pm 0,1; 50-100 Kg \pm 0,2;
100-160 Kg \pm 0,3

This scale cannot be used for certified weight measurements. We reserve the right to apply technical changes.

WARRANTY

This device is covered by a 2-year warranty from the date of purchase, against manufacturing defects, valid only upon submission of the purchase receipt. Opening or tampering with the device voids the warranty. The warranty does not cover accidental damage caused by improper use or failure to abide by the instructions for use.

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