**COLPHARMA**  <sup>®</sup>

*Caldo Piu'*

Heating pad

**12173**  
**BASIC model**

**12174**  
**BACK & NECK model**


EN

IT

DE

FR

**Temperature selector / Selettore di temperatura / Temperaturwahlschalter /  
Réglage température**

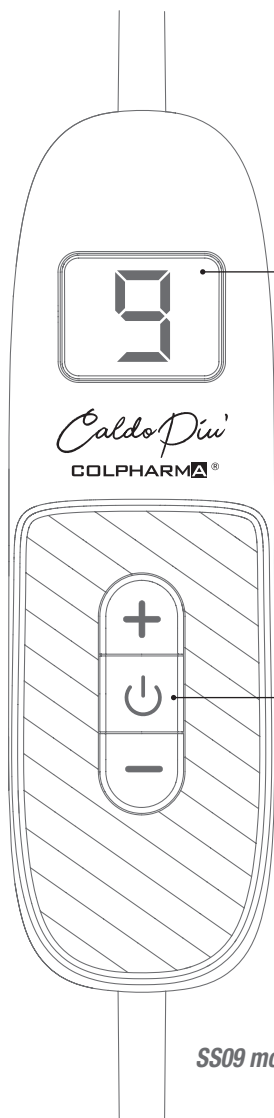9-level temperature display  
Display 9 livelli di temperatura  
Display mit 9 Temperaturstufen  
Affichage 9 niveaux de température

ON / OFF button  
Tasto ON / OFF  
EIN / AUS - Taste  
Bouton ON/OFF

**SS09 model**

# Easy to clean! / Facile da pulire! / Leicht zu reinigen! / Facile à nettoyer !

## Washing instructions / Istruzioni per il lavaggio / Waschanleitung / Instructions pour le lavage

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## User's Instructions

**Read this manual carefully before use.**

**Please retain the user's instructions for future reference.**

The CaldoPiù heating pad offers convenience and ease of use, thanks to its flexible and comfortable structure, as well as its washing machine compatibility.

It heats up quickly, ensuring even and consistent warmth. The 9 temperature levels allow effortless adjustment of the heat intensity throughout the therapy session. It is safe to use and includes an auto shut-off feature that activates after 90 minutes of continuous use.

Utilizing the heating pad provides an immediate sense of well-being, alongside its analgesic and pain-relieving effects.

## WARNING

Please carefully review the user's instructions, specifically those pertaining to safety, before operating the device. It is essential to retain these instructions for future reference.

In the event that the device is transferred to another person, it is imperative to accompany it with these user's instructions. This heating pad is specifically designed to emit therapeutic heat and adheres to international standards for such products.

The perception of heat emitted by the heating pad may vary when it comes into contact with the skin. While the pad may feel lukewarm when touched in an open environment, the sensation of heat will intensify when worn.

To achieve your optimal temperature, adjust the heating pad accordingly while wearing it.

## USE PRECAUTIONS

- This device must be used exclusively as described in this manual. The manufacturer/ distributor shall not be held liable for any damages resulting from improper use.
- Before each use, carefully inspect the heating pad. Do not use it if you notice any signs of wear or damage.
- Prior to plugging the device into a power outlet, ensure that the nominal voltage on the nameplate is compatible with the electrical grid in your location.
- This device contains delicate components and should be handled with care. Follow the storage and cleaning guidelines outlined in the "Technical Specifications" and "Device Cleaning" sections.
- Do not use the device on children, disabled individuals, sleeping individuals, individuals with reduced

- sensitivity, unconscious individuals, diabetic individuals, or individuals with poor blood circulation.
- The device can be used by children aged 8 and above under adult supervision or after receiving instructions on safe usage and understanding the associated risks.
  - Use the device while lying down, placing it on the desired body part for warming.
  - Prolonged use of the heating pad at high temperature settings may cause skin irritation or burns.
  - For prolonged use, it is recommended that the temperature be set to “1”.
  - Do not use the heating pad on swollen or injured body parts.
  - Do not fall asleep while the heating pad is in use.
  - Do not use the device if you are in contact with water or other liquids. Only use it in dry environments (avoiding bathrooms or similar locations).
  - In the event of the device accidentally falling into water, immediately unplug it from the power outlet and contact customer service before reuse.
  - Do not puncture the heating pad or use pins, needles, or any other sharp objects to position the device as desired.
  - Avoid folding the heating pad onto itself while it is operational.
  - Do not sit or lie on the device and avoid crushing or crumpling it.
  - Do not cover the heating pad with pillows, throws, blankets, or any other items. The temperature selector should not be covered, placed above, or placed below the heating pad while the device is in operation.
  - Do not use the device as a bed warmer or blanket warmer during sleep.
  - Keep the device and power cord away from heat sources.
  - Never transport, pull, or twist the heating pad by the power cord, and ensure that the cord does not become tangled.
  - Avoid applying creams or ointments that may increase heat.
  - Regularly check the skin in contact with the device to prevent irritation or burns.
  - Always disconnect the device from the power outlet after use, without exception.
  - Do not hand wash or machine wash the device unless the power plug is disconnected from the power outlet and the heating pad is detached from the power cord connector.
  - Do not use a tumble dryer to dry the heating pad, and do not iron it.
  - The heating pad can only be used with the provided temperature selector (SS09).
  - Allow the heating pad to cool before folding it for storage. The manufacturer accepts liability for the safety, reliability, and performance of the device only when the electrical installation of the premises where the device is used complies with legal regulations and the device is used as recommended in this user’s instructions.
  - This device is not suitable for hospital use or in oxygen-enriched environments.
  - If you experience pain or discomfort while using the device, immediately turn off the heating pad.
  - Consult a physician before using the device if you have any health concerns.
  - If persistent muscle and joint pain occurs, inform your physician.

## USING THE DEVICE FOR THE FIRST TIME

---

- Insert the device plug into a 230V electrical outlet.
- Power on the device by pressing the ON/OFF «  » button located on the temperature selector.
- Adjust the temperature selector to the desired level (ranging from 1 to 9) using the “+” and “-” buttons.

- Once the desired heat is achieved, apply the device to the intended body area.
- When finished, turn off the device by pressing the ON/OFF «  » button on the temperature selector.
- Disconnect the plug from the electrical outlet.
- Allow the device to cool down and store it in its original packaging.

## DEVICE CLEANING

---

- Prior to cleaning the device, unplug it from the electrical outlet and detach the connector from the heating pad. Let it cool for a minimum of 10 minutes.
- Wash the device in a washing machine at a maximum temperature of 30°C or clean it by hand using a mild detergent. Alternatively, you can wipe the device with a damp cloth without using any detergents or cleaning agents.
- Place the device in a clean and dry area, laying it flat without any covering or items resting on top. Ensure that it is completely dry before reattaching the connector to the heating pad and plugging it back into the electrical outlet.

**The device must be disposed of according to local regulations and not as ordinary household waste.**

## TECHNICAL SPECIFICATIONS

---

<b>Model 12173 dimensions:</b>	40x30 cm
<b>Model 12174 dimensions:</b>	42x63 cm
<b>Voltage:</b>	220-240V~50/60 Hz
<b>Power:</b>	100W
<b>Heating levels:</b>	from 1 to 9 levels
<b>Automatic shut-off:</b>	after 90 minutes
<b>Max. temperature:</b>	70±5 °C

We reserve the right to apply technical changes

## WARRANTY

---

This device is covered by a **2-year** warranty from the date of purchase, against manufacturing defects, valid only upon submission of the purchase receipt. Opening or tampering with the device voids the warranty. The warranty does not cover accidental damages caused by improper use or failure to follow the user instructions.

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